

DTIC FILE COPY

2

TECHNICAL REPORT BRL-TR-3163

BRL

AD-A228 373

LASER MICROPLASMA-GAS CHROMATOGRAPHY DETECTOR.
II. DETECTION OF SPECIES-SPECIFIC FRAGMENT EMISSION

RANDY J. LOCKE
JEFFREY B. MORRIS
BRAD E. FORCH
ANDRZEJ W. MIZIOLEK

OCTOBER 1990

DTIC
ELECTE
NOV 01 1990
S B D

APPROVED FOR PUBLIC RELEASE; DISTRIBUTION UNLIMITED.

U.S. ARMY LABORATORY COMMAND

BALLISTIC RESEARCH LABORATORY
ABERDEEN PROVING GROUND, MARYLAND

90 10 21 019

NOTICES

Destroy this report when it is no longer needed. DO NOT return it to the originator.

Additional copies of this report may be obtained from the National Technical Information Service, U.S. Department of Commerce, 5285 Port Royal Road, Springfield, VA 22161.

The findings of this report are not to be construed as an official Department of the Army position, unless so designated by other authorized documents.

The use of trade names or manufacturers' names in this report does not constitute indorsement of any commercial product.

UNCLASSIFIED

REPORT DOCUMENTATION PAGE			Form Approved OMB No. 0704-0188
<small>Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503.</small>			
1. AGENCY USE ONLY (Leave blank)	2. REPORT DATE October 1990	3. REPORT TYPE AND DATES COVERED Final, Nov 89 - Aug 90	
4. TITLE AND SUBTITLE LASER MICROPLASMA-GAS CHROMATOGRAPHY DETECTOR. II. DETECTION OF SPECIES-SPECIFIC FRAGMENT EMISSION		5. FUNDING NUMBERS 1L161102AH43	
6. AUTHOR(S) Randy J. Locke, Jeffrey B. Morris, Brad E. Forch, Andrzej W. Miziolek			
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES)		8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES) Ballistic Research Laboratory ATTN: SLCBR-DD-T Aberdeen Proving Ground, MD 21005-5066		10. SPONSORING/MONITORING AGENCY REPORT NUMBER BRL-TR-3163	
11. SUPPLEMENTARY NOTES To be published in 22 November 1990 edition of <u>Applied Optics</u> .			
12a. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release; distribution unlimited		12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words) Characteristic laser-produced microplasma emissions from various simple carbon-containing vapors entrained in a He carrier gas have been observed and compared. A focused ArF (193 nm) excimer laser is used to induce microplasmas with modest pulse energies (15 mJ or less) in the effluent region of a gas chromatography (GC) capillary column. Strong atomic (C, H, O, Cl, and F), as well as molecular (C ₂ ⁺ , CH, and CCl) emissions are observed with very high signal-to-noise ratios. A plasma emission survey indicates that different classes of molecules show unique spectra which make it relatively easy to distinguish one chemical class from another. These results suggest that a laser microplasma gas chromatography detector (LM-GCD) should offer additional discrimination/resolution for unknown sample gas mixture analysis. In addition, the LM-GCD exhibits a significant advantage over certain other GC detectors, like the widely used flame ionization detector (FID), by readily detecting non-responsive gases such as CO, CO ₂ , CCl ₄ , and freons. (72) ←			
14. SUBJECT TERMS Laser-Produced Microplasmas, Gas Chromatography Detector, Plasma Emission; Gas Chromatography		15. NUMBER OF PAGES 26	
		16. PRICE CODE	
17. SECURITY CLASSIFICATION OF REPORT Unclassified	18. SECURITY CLASSIFICATION OF THIS PAGE Unclassified	19. SECURITY CLASSIFICATION OF ABSTRACT Unclassified	20. LIMITATION OF ABSTRACT UL

INTENTIONALLY LEFT BLANK.

TABLE OF CONTENTS

	<u>Page</u>
ACKNOWLEDGEMENTS	v
I. INTRODUCTION	1
II. EXPERIMENTAL	1
III. RESULTS AND DISCUSSION	4
IV. CONCLUSION	10
REFERENCES	13
DISTRIBUTION LIST	15



Accession For	
NTIS GRA&I	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification _____	
By _____	
Distribution/	
Availability Codes	
Dist.	Avail and/or Special
A-1	

INTENTIONALLY LEFT BLANK

ACKNOWLEDGEMENTS

We would like to thank Michael Ellzy and Monica Heyl at CRDEC for supplying a Megabore column. The authors are also grateful to Dr. Stephen Howard for his assistance in manuscript preparation. RJI would like to acknowledge the BRL/NRC postdoctoral program for support. This work was supported, in part, by the Air Force Office of Scientific Research, Directorate of Aerospace Sciences, Contract # 90-0025.

INTENTIONALLY LEFT BLANK

I. INTRODUCTION

Pulsed lasers have been used to produce microplasmas via laser-induced breakdown almost from the time of their invention. Virtually all of the reports in the intervening years involve the use of the ruby, Nd:Glass/Nd:YAG, or CO₂ lasers.¹⁻³ More recently, a specific technique called LIBS (laser-induced breakdown spectrometry) was developed for the purpose of chemical (elemental) analysis of aerosol and gas samples.^{4,5} In the last few years a new approach for laser plasma formation in gases has been developed which utilizes ultraviolet lasers and is based on atomic and molecular resonances in the microplasma formation processes.^{6,7} Attributes of this approach include efficiency with respect to the amount of incident laser energy that is required to form microplasmas, i.e., a lower threshold for breakdown, as well as a much higher degree of control over the intensity of the microplasmas, especially near the threshold region. This latter effect allows for the formation of much weaker microplasmas using UV resonant lasers than that possible with visible/infrared non-resonant lasers. As such, the weaker (lower temperature) microplasmas contain a relatively larger fraction of molecular fragments and, conversely, a smaller fraction of atoms and ions.

We recently reported our first results of the use of resonant laser-produced microplasmas for chemical analysis in a study which was aimed at ascertaining the utility of the laser microplasma approach for potential use as a detector for gas chromatography.⁸ Very simply, an ArF excimer laser beam (193 nm) was focused at the orifice of a gas chromatography (GC) capillary column through which a high purity He carrier gas flowed continuously. Utilizing modest laser pulse energies (ca. 10 mJ), microplasmas were formed when a small amount of analyte sample was injected into the carrier gas flow and were not produced in the He carrier gas by itself. Resonance in the formation of free electrons in the laser focal volume was accomplished by the coincident absorption of photons from the broadband excimer laser by carbon atoms at 193.1 nm ($2p^2\ ^1D_2 \rightarrow 2p3s\ ^1P_1^o$).⁸ Microplasma formation was detected using photoacoustic, optogalvanic, and photometric emission spectroscopy. Photometric detection was found to be the most sensitive method of the three techniques and thus offered the greatest promise for further development.

In this paper we report further work on the development of this technique which involves a survey of spectrally-resolved microplasma emissions for a variety of carbon-containing sample gases which are meant to represent different types of chemical groups. The results indicate that there is significant group-specificity in the emission spectra which should make it possible to increase the chemical selectivity/resolution of a gas chromatograph based on microplasma emission analysis using the LM-GCD technique.

II. EXPERIMENTAL

The apparatus and techniques used for the study of the dispersed emission of electronically excited fragments formed within the microplasma of various gaseous flows are essentially identical to those previously employed by Morris, et al.⁸ A schematic illustration of the experimental arrangement is presented in Figure 1A. The excitation

source is a Lumonics Model TE 861M-3 pulsed, broadband ArF (193 nm) excimer laser, which is equipped with unstable resonator optics. The laser pulse duration is approximately 15 nsec and, for this study, pulse energies of up to 15 mJ are used. The laser is operated at a repetition rate of 15 Hz. The laser pulses are focused through a 25 mm diameter, 50 mm focal length UV grade fused silica lens. The focus of the beam is directed to a position approximately 0.5 mm above the outlet of a 5 meter, 544 micron diameter fused silica "Megabore" column obtained from J & W Scientific which is normally used in the HP Model 5890A gas chromatograph. Since the purpose of this present study is to characterize microplasma emission from individual analytes and not analyte mixtures, the GC column was removed from the GC and used separately for experimental convenience. The laser energy is monitored by a Scientech laser power meter which also serves as a beam stop.

The LM-GCD comprises two separate components, specifically the sheath gas and column effluent. To prevent the mixing of room air with the GC effluent gases, the column outlet, illustrated in Figure 1B, is encircled by a 3/8" diameter helium sheath with a flow rate of approximately 200 ml/minute. The use of a sheath gas is unique to the LM-GCD. One important reason for separating room air and the GC effluent gases is to prevent a condition of breakdown in 100% of the laser shots. This is particularly important when using the LM-GCD for maximum sensitivity where broadband, rather than wavelength-specific, detection of weak microplasma emission signifies the elution of a small quantity of the analyte molecules.⁸ The room air-induced breakdown is due to the ease with which molecular oxygen is ionized by the ArF excimer laser via the resonant two-photon absorption through the intermediate Schumann-Runge (S-R) bands.⁹ Although absorption by laboratory air will clearly deplete the laser radiation in the O₂ S-R band wavelength region, at short distances from the laser (approximately 0.5 meter in our case), this effect is minimal.⁹ Another reason for excluding room air from the laser focal volume is to prevent the contamination of the microplasma emission by gases which are not components of the original sample. Centering of the column outlet within the sheath gas is accomplished by the utilization of a wire mesh which also allows for the unrestricted flow of the sheath gas.

The column effluent gas is composed of carrier gas and analyte vapor. The helium carrier gas (Spectra Gases 99.999%), after successive filtering and drying, flows to the atmosphere through the Megabore column at a rate of roughly 20 ml/minute. Gaseous samples were introduced into the carrier gas prior to the Megabore column at a typical flow rate of also about 20 ml/minute resulting in approximately 10¹³ analyte molecules present in the laser probe volume (10⁻⁶ cm³). Backing pressures and flow controls were adjusted to maintain similar concentrations in the effluent of the capillary column for the different sample gases. Alternately, liquid phase compounds (for most samples the vapor pressures at 298 K are in the 100 torr range) were inserted into the column by flowing the carrier gas over the neat liquid. This resulted in small variations of sample vapor concentrations in the probe volume between the different liquids. The difficulty of maintaining similar concentrations for the different liquid samples precluded the quantification and absolute comparison of the emission results.

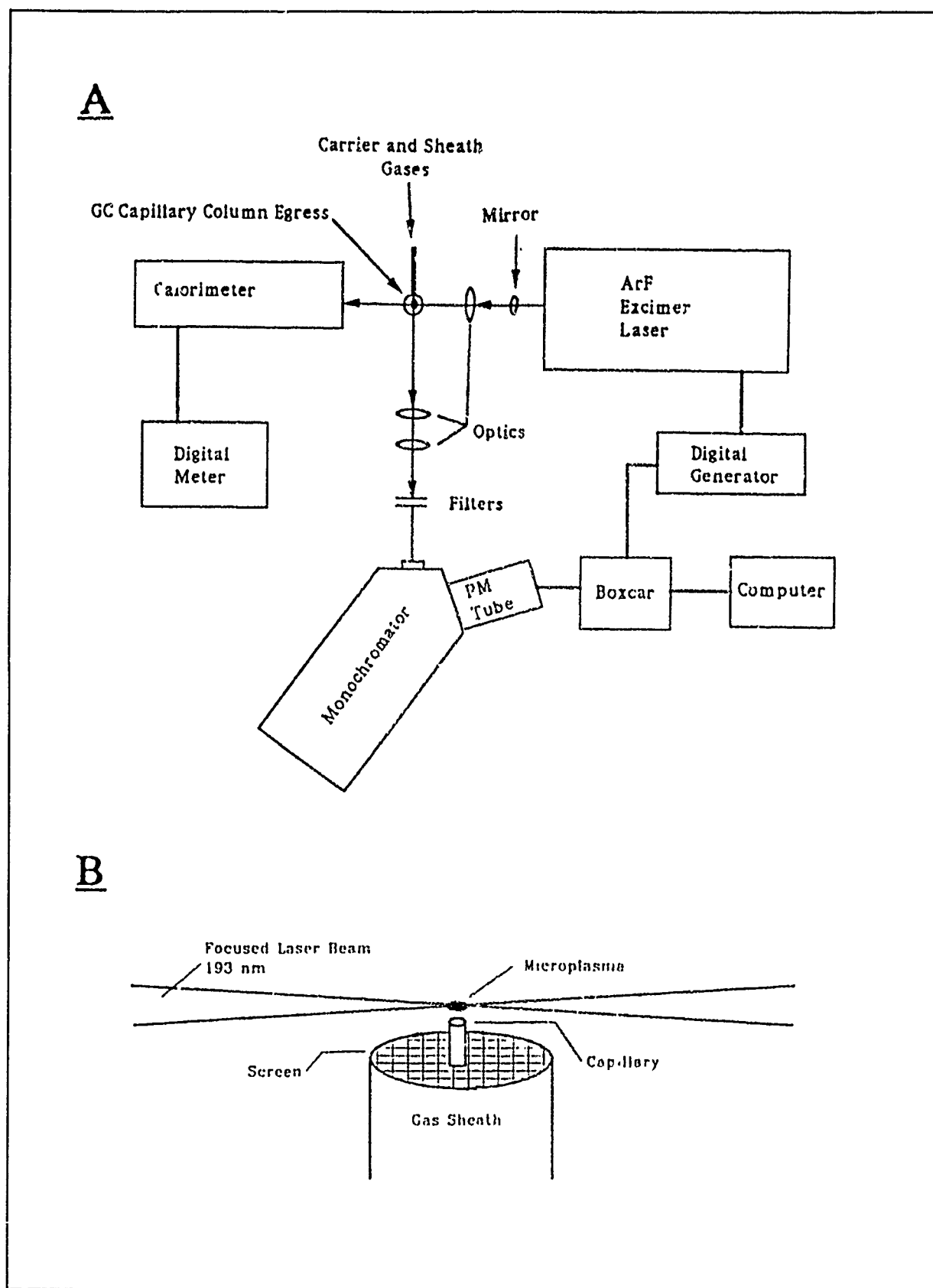


Figure 1. (A) Schematic of Experimental Apparatus. (B) Schematic Enlargement of Laser Interaction Region with the Effluent Gases from the Capillary Column.

A two-lens collection system, perpendicular to the direction of the laser excitation source, focused the emission originating within the microplasmas onto the entrance slit (width = 250 microns) of a 0.25 meter McPherson monochromator equipped with a Hamamatsu 955R photomultiplier tube (PMT). The signal from the PMT was processed with a Stanford Research Systems boxcar integrator with a gate width of 100 nsec. The output of the boxcar was then recorded utilizing a Zenith 80286-based computer. The spectra were acquired in three segments of approximately equal wavelength span. This was necessary in order to minimize the drop in ArF laser energy, and thus the signal level, over a prolonged spectral scan. This behavior, unfortunately, is quite characteristic of older excimer lasers such as ours. Typically these scans encompassed the following regions; 190-400 nm, 400-600 nm, and 600-850 nm. Consolidation of these three spectral regions into one spectrum was accomplished through use of the 1-2-3 Access System from Lotus Development Corporation. It should be noted that prior to acquiring the longer wavelength scans (>400 nm), a filter (GG-400) was inserted into the dispersed emission beam path in order to eliminate appearance of short wavelength emission in second order of grating from being recorded.

III. RESULTS AND DISCUSSION

To demonstrate the versatility and general nature of the LM-GCD, a number of different carbon containing compounds were examined. Figures 2, 3 and 4 present representative emission spectra originating from the electronically excited fragments formed within the various microplasmas. Additionally, Table 1 provides a detailed assignment of the major peaks observed in each respective spectrum. The relative intensity values are uncorrected for the monochromator/PMT instrument response function and the wavelengths listed represent actual readings from the monochromator mechanical dial (accuracy ± 0.5 nm).

Figure 2 displays the emission spectra obtained from the microplasmas generated in a flow of carbon monoxide (CO) and methanol (CH₃OH), top and bottom, respectively. Both CO and CH₃OH exhibit strong atomic carbon and oxygen emissions. CH₃OH, however, displays the additional emissions attributed to the atomic H α and H β transitions. It is noteworthy to consider that had a conventional flame ionization detector (FID) been utilized, only the presence of CH₃OH would have been discernible. The non-combustible CO would have remained undetected. The carbon dioxide (CO₂) microplasma emission spectrum was found to be very similar to that of CO (see Table 1). It should be noted that all of the spectra acquired for Table 1 display high signal-to-noise ratios for most atomic and molecular emissions. In addition, it was determined that the broad band appearing at 387 nm in this and all subsequent spectra is attributed to second order ArF laser scattering.

Figure 3 shows the emission spectra acquired from electronically excited fragments produced in a microplasma generated in a flow of propane (CH₃CH₂CH₃) and benzene (C₆H₆), top and bottom, respectively. It is immediately apparent that the only feature common to both compounds is that ascribed to the emission of atomic carbon (448 nm). Two additional lines observed in the propane spectrum are those of H α and H β .

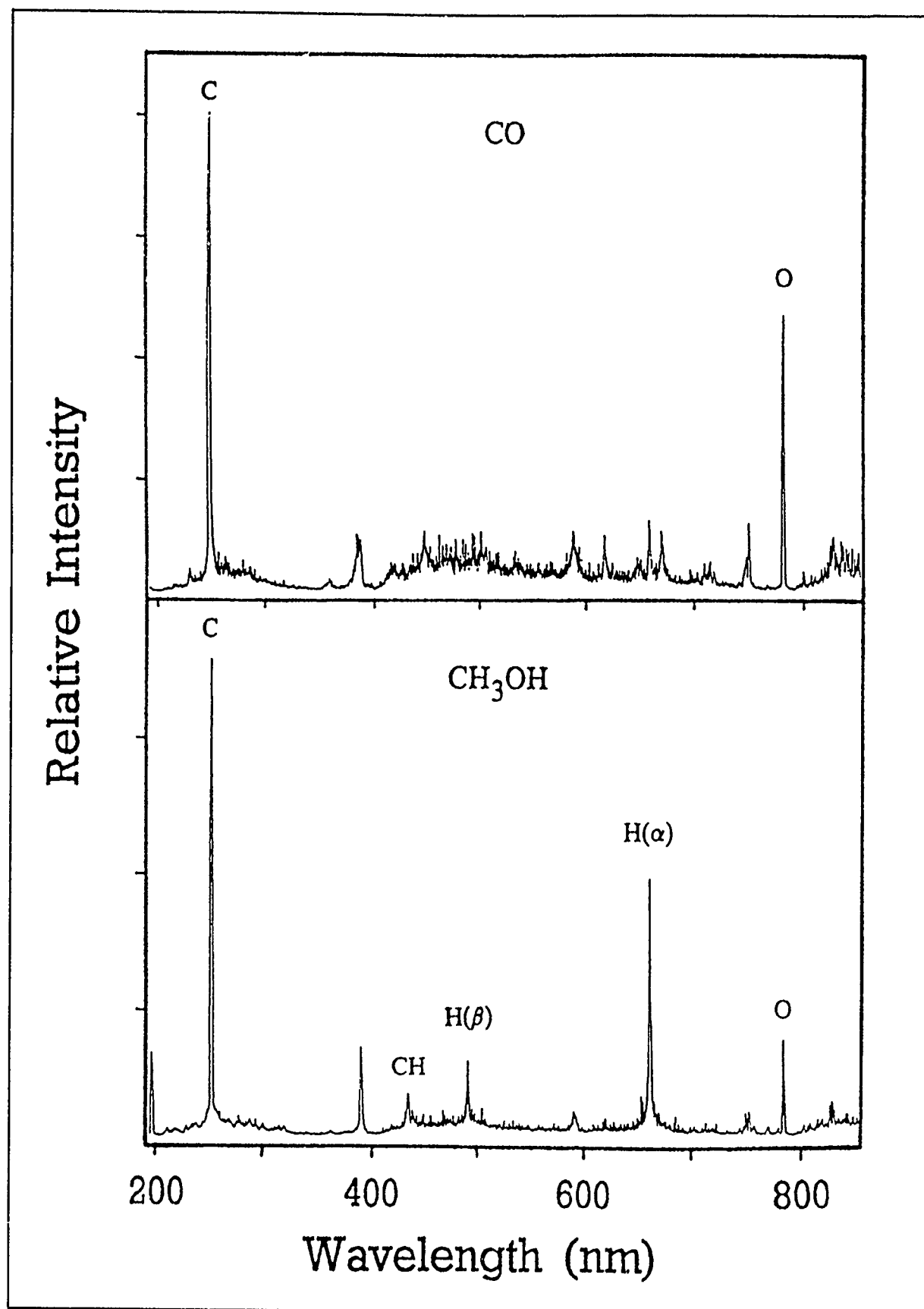


Figure 2. ArF Laser-Produced Microplasma Emission Spectra for CO and CH₃OH

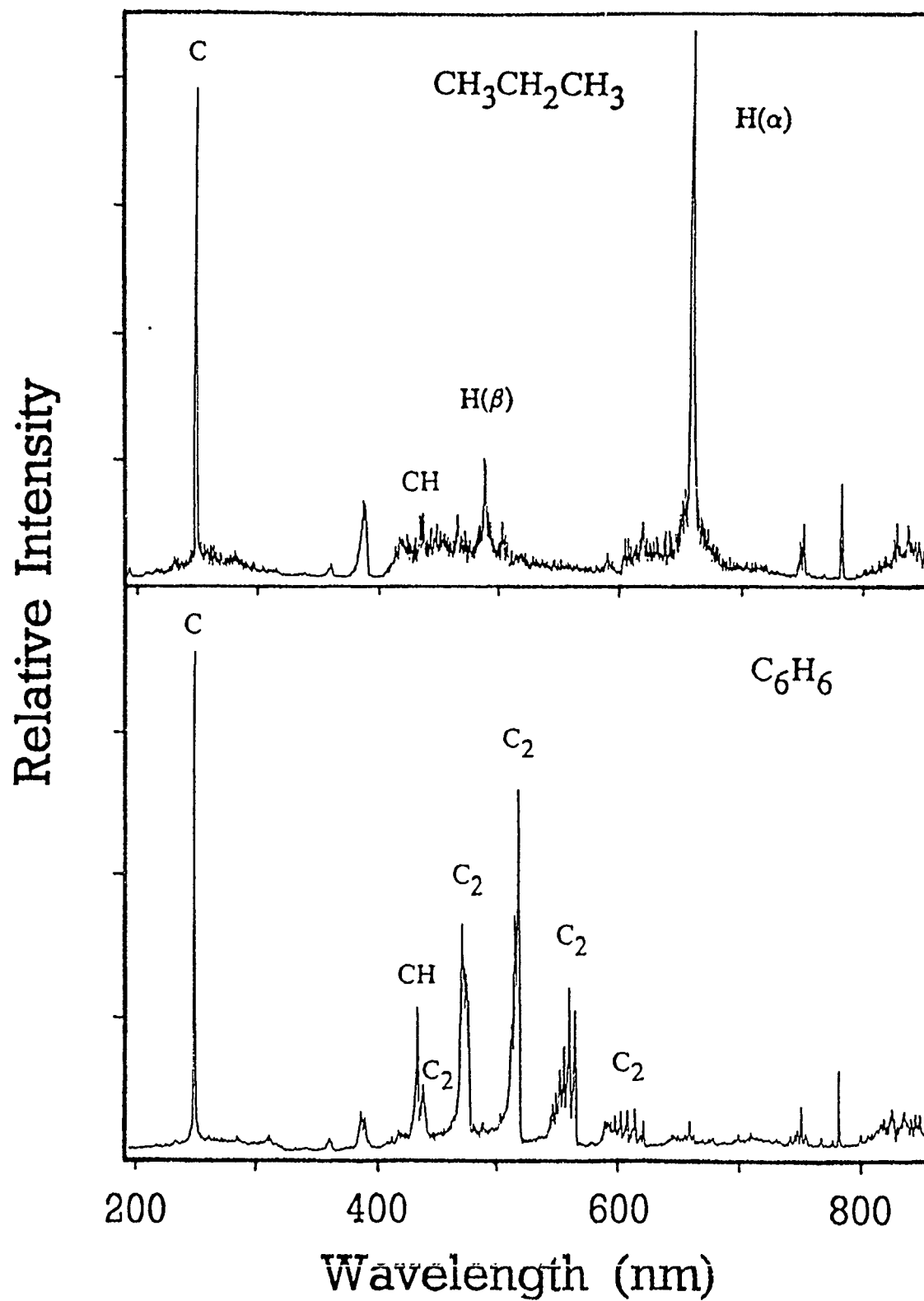


Figure 3. ArF Laser-Produced Microplasma Emission Spectra for Propane and Benzene

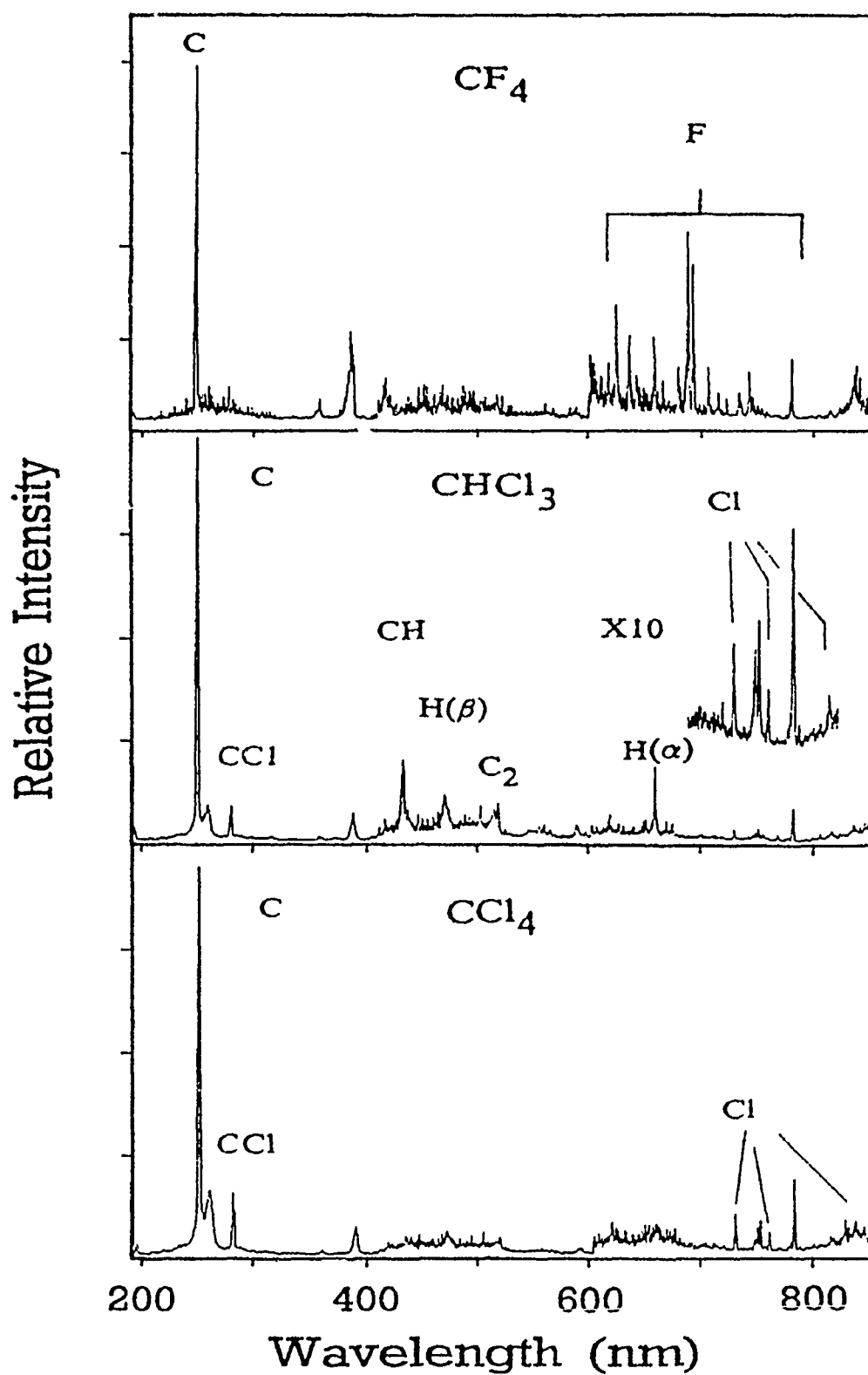


Figure 4. ArF Laser-Produced Microplasma Emission Spectra for CF_4 , CHCl_3 , and CCl_4

Table 1. Laser-Produced Microplasma Emissions from Various Carbon Containing Sample Molecules

Molecule	Overall Plasma I ^a	Species ^b	Emission nm ^c	I ^d
CH ₃ CH ₂ CH ₃	1.00	C	248.2	7.8
		CH	431.8	1.1
		H(β)	486.1	2.0
		H(α)	657.2	8.8
		O	778.2	2.6
CH ₄	0.74	C	248.3	5.8
		CH	432.0	1.2
		H(β)	486.5	3.8
		C ₂ (0,0)	516.7	0.9
		H(α)	657.2	6.7
CH ₃ OH	e	C	248.2	8.0
		CH	431.3	0.8
		H(β)	487.0	1.3
		H(α)	657.0	4.3
		O	777.8	1.5
CO	0.92	C	248.3	7.2
		O	778.2	4.3
CO ₂	0.56	C	248.3	4.4
		O	777.8	2.5
C ₆ H ₆	e	C	248.3	7.1
		CH	431.5	2.1
		C ₂ (2,0)	438.1	1.0
		C ₂ (1,0)	474.1	2.2
		C ₂ (0,0)	516.9	5.1
		C ₂ (0,1)	564.0	2.3
		C ₂ (0,2)	619.6	0.5
		H(α)	657.4	0.6
		O	778.2	1.3
C ₆ H ₅ F	e	C	248.2	10.0
		CH	431.2	4.0
		C ₂ (2,0)	437.8	2.6
		C ₂ (1,0)	474.0	5.2
		C ₂ (0,0)	516.5	9.2
		C ₂ (0,1)	563.8	4.6
		C ₂ (0,2)	619.2	1.4
		H(α)	657.0	1.1

Table 1. Laser-Produced Microplasma Emissions from Various Carbon Containing Sample Molecules				
Molecule	Overall Plasma I ^a	Species ^b	Emission nm ^c	I ^d
CHCl ₃	e	C	248.3	7.9
		CCl	278.3	0.8
		CH	431.3	1.7
		H(β)	468.5	0.9
		C ₂ (0,0)	516.7	0.8
		H(α)	657.0	1.5
		Cl	726.4	0.3
		O	755.3	0.2
			809.2	0.1
			778.0	0.7
Cl ₄	e	C	248.3	9.6
		CCl	278.3	1.6
		C ₂ (0,0)	516.8	0.5
		Cl	725.8	1.1
			755.5	0.5
			809.6	0.4
			778.2	2.0
CF ₄	0.67	C	248.3	5.2
		F	623.8	1.3
			634.9	0.8
			641.5	0.9
			686.1	2.8
			691.1	2.2
			704.8	0.9
			740.7	0.8
		O	777.7	1.0

Annotations:

- Plasma emission intensities relative to strongest carbon emission (CH₃CH₂CH₃).
- The transitions listed for C₂ are (ν' , ν'').
- Approximate wavelength readings obtained directly from monochromator.
- Intensity is given in arbitrary units and is uncorrected for instrument response.
- Liquid samples (see text).

Conversely, the spectrum of benzene is dominated by the presence of a complex series of band structures. These features are found to be the various transitions of the C_2 Swan system of bands ($A^3\Pi_g \leftarrow X^3\Pi_u$);¹⁰ (2,0) at 438.1 nm, (1,0) at 474.1 nm, (0,0) at 516.9 nm, (0,1) at 564.0 nm, and (0,2) at 619.6 nm. Additional emissions noted in the benzene microplasma spectrum are CH and H_α . Spectral results similar to those obtained for benzene and propane were also obtained from microplasmas formed in the respective flows of fluorobenzene (C_6H_5F) and methane (CH_4) (see Table 1). The emission observed at approximately 777.8 nm in these spectra is attributed to atomic oxygen and is indicative of the difficulty of totally isolating the GC effluent gases from room air. One solution to this problem under evaluation is the construction of a flow cell in which the effluent gases can be probed without interference from room air.

Figure 4 displays the emission spectra of electronically excited fragments from the microplasmas generated in the flows of Freon-14 (CF_4), chloroform ($CHCl_3$) and carbon tetrachloride (CCl_4), top, middle and bottom, respectively. Upon close inspection of Figure 4, the only common feature is found to be that attributed to the emission of atomic carbon. A closer examination of the $CHCl_3$ spectrum reveals that in addition to the CH bands and the H_α and H_β lines, a new set of lines appears in the 700-825 nm region. These lines, which are assigned to the emission of atomic chlorine,^{5,11} also appear in the microplasma spectrum of CCl_4 . Furthermore, a band visible at 278.3 nm in the emission spectra of only $CHCl_3$ and CCl_4 is attributed to the $CCl A^2\Delta \leftarrow X^2\Pi$ (0,0) transitions.¹⁰ The identity of the other band at 257.7 nm common to both $CHCl_3$ and CCl_4 is as of yet unresolved.

The longer wavelength region of CF_4 (600-800 nm) is remarkable due to the presence of an elaborate set of relatively intense lines wholly attributable to atomic fluorine emission.^{5,11} As with the previous examples, Figure 4 demonstrates the species specificity of the laser-produced microplasma emission technique. Of the three compounds shown in Figure 4, conventional FID is unresponsive to both CF_4 and CCl_4 .

IV. CONCLUSION

It has long been known that many substances are not amenable to detection by conventional GC methods. Those compounds that produce little or no response in FID (e.g., O_2 , H_2 , CO, CO_2 , halogenated compounds, etc.) must therefore have some other means for detection. Presented here are qualitative laser-produced microplasma emission results which indicate that the LM-GCD offers not only broad chemical response, but also considerable chemical specificity. These results, in addition to our previous work which demonstrated considerable sensitivity potential, suggest that the LM-GCD approach is worthy of serious consideration as a universal carbon detector.

Further improvements are currently being made to the experimental apparatus. These include the replacement of the photomultiplier with an optical multichannel analyzer (OMA) which will allow for quicker data collection and the ability to take large portions of the spectrum simultaneously. Future work will concentrate on the time resolution aspects of the microplasma emission as well as on the determination of the realistic

detection limits of which the LM-GCD is capable. In general, a number of technological advances should make it possible for practical implementation of the LM-GCD. These include the development of a compact and inexpensive UV laser* as well as improvement in the UV transmission characteristics of optical fibers.

*A contract to design and construct a compact and field-worthy UV laser under the BRL SBIR Phase I program is in progress.

INTENTIONALLY LEFT BLANK

REFERENCES

1. C.G. Morgan, "Laser-Induced Breakdown of Gases", Rep. Prog. Phys., Vol. 38, p. 621, 1975.
2. D.C. Smith and R.G. Meyerand, Jr., "Laser Radiation Induced Gas Breakdown", Princ. Laser Plasmas, G. Bekefi, ed., p. 457, Wiley, New York, 1976.
3. P.K. Carroll and E.T. Kennedy, "Laser-Produced Plasmas", Contemp. Phys., Vol. 22, p. 61, 1981.
4. L.J. Radziemski, T.R. Loree, D.A. Cremers, and N.M. Hoffman, "Time-Resolved Laser-Induced Breakdown Spectrometry of Aerosols", Anal. Chem., Vol. 55, p. 1246, 1983.
5. D.A. Cremers and L.J. Radziemski, "Detection of Chlorine and Fluorine in Air by Laser-Induced Breakdown Spectrometry", Anal. Chem., Vol. 55, p. 1252, 1983.
6. B.E. Forch and A.W. Miziolek, "Oxygen-Atom Two-Photon Resonance Effects in Multiphoton Photochemical Ignition of Premixed H_2/O_2 Flows", Opt. Lett., Vol. 11, p. 129, 1986.
7. B.E. Forch and A.W. Miziolek, "Ultraviolet Laser Ignition of Premixed Gases by Efficient and Resonant Multiphoton Photochemical Formation of Microplasmas", Combust. Sci. & Tech., Vol. 52, p. 151, 1987.
8. J.B. Morris, B.E. Forch, and A.W. Miziolek, "A Novel Detector for Gas Chromatography Based on UV Laser-Produced Microplasmas", Appl. Spectrosc., Vol. 44, No. 6, p. 1040, 1990.
9. M.P. Lee and R.K. Hanson, "Calculations of O_2 Absorption and Fluorescence at Elevated Temperatures for a Broadband Argon-Fluoride Laser Source at 193 nm", J. Quant. Spectrosc. Radiat. Transfer, Vol. 36, p. 425, 1986.
10. R.W. Pearce and A.G. Gaydon, The Identification of Molecular Spectra, 4th ed., Halsted Press, London, 1976.
11. Atomic Transition Probabilities, Vols. 1 and 2, National Standard Reference Data Series 4, National Bureau of Standards, 1966.

INTENTIONALLY LEFT BLANK

<u>No. of Copies</u>	<u>Organization</u>
1	Office of the Secretary of Defense OUSD(A) Director, Live Fire Testing ATTN: James F. O'Bryon Washington, DC 20301-3110
2	Administrator Defense Technical Info Center ATTN: DTIC-DDA Cameron Station Alexandria, VA 22304-6145
1	HQDA (SARD-TR) WASH DC 20310-0001
1	Commander US Army Materiel Command ATTN: AMCDRA-ST 5001 Eisenhower Avenue Alexandria, VA 22333-0001
1	Commander US Army Laboratory Command ATTN: AMSLC-DL Adelphi, MD 20783-1145
2	Commander US Army, ARDEC ATTN: SMCAR-IMI-I Picatinny Arsenal, NJ 07806-5000
2	Commander US Army, ARDEC ATTN: SMCAR-TDC Picatinny Arsenal, NJ 07806-5000
1	Director Benet Weapons Laboratory US Army, ARDEC ATTN: SMCAR-CCB-TL Watervliet, NY 12189-4050
1	Commander US Army Armament, Munitions and Chemical Command ATTN: SMCAR-ESP-L Rock Island, IL 61299-5000
1	Commander US Army Aviation Systems Command ATTN: AMSAV-DACL 4300 Goodfellow Blvd. St. Louis, MO 63120-1798

<u>No. of Copies</u>	<u>Organization</u>
1	Director US Army Aviation Research and Technology Activity
1	ATTN: SAVRT-R (Library) M/S 219-3 Ames Research Center Moffett Field, CA 94035-1000
1	Commander US Army Missile Command ATTN: AMSMI-RD-CS-R (DOC) Redstone Arsenal, AL 35898-5010
1	Commander US Army Tank-Automotive Command ATTN: AMSTA-TSL (Technical Library) Warren, MI 48397-5000
1	Director US Army TRADOC Analysis Command ATTN: ATAA-SL White Sands Missile Range, NM 88002-5502
(Class. only) 1	Commandant US Army Infantry School ATTN: ATSH-CD (Security Mgr.) Fort Benning, GA 31905-5660
(Unclass. only) 1	Commandant US Army Infantry School ATTN: ATSH-CD-CSO-OR Fort Benning, GA 31905-5660
1	Air Force Armament Laboratory ATTN: AFATL/DLODL Eglin AFB, FL 32542-5000
	<u>Aberdeen Proving Ground</u>
2	Dir, USAMSAA ATTN: AMXSY-D AMXSY-MP, H. Cohen
1	Cdr, USAFTECOM ATTN: AMSTE-TO-F
3	Cdr, CRDEC, AMCCOM ATTN: SMCCR-RSP-A SMCCR-MU SMCCR-MSI
1	Dir, VLAMO ATTN: AMSLC-VL-D

No. of Copies	Organization
4	Commander US Army Research Office ATTN: R. Ghirardelli D. Mann R. Singleton R. Shaw P.O. Box 12211 Research Triangle Park, NC 27709-2211
2	Commander US Army, ARDEC ATTN: SMCAR-AEE-B, D.S. Downs SMCAR-AEE, J.A. Lannon Picatinny Arsenal, NJ 07806-5000
1	Commander US Army, ARDEC ATTN: SMCAR-AEE-BR, L. Harris Picatinny Arsenal, NJ 07806-5000
2	Commander US Army Missile Command ATTN: AMSMI-RK, D.J. Ishin W. Wharton Redstone Arsenal, AL 35898
1	Commander US Army Missile Command ATTN: AMSMI-RKA, A.R. Maykut Redstone Arsenal, AL 35898-5249
1	Office of Naval Research Department of the Navy ATTN: R.S. Miller, Code 432 800 N. Quincy Street Arlington, VA 22217
1	Commander Naval Air Systems Command ATTN: J. Ramnarace, AIR-54111C Washington, DC 20360
1	Commander Naval Surface Warfare Center ATTN: J.L. East, Jr., G-23 Dahlgren, VA 22448-5000
2	Commander Naval Surface Warfare Center ATTN: R. Bernecker, R-13 G.B. Wilnot, R-16 Silver Spring, MD 20903-5000

No. of Copies	Organization
5	Commander Naval Research Laboratory ATTN: M.C. Lin J. McDonald E. Oran J. Shnur R.J. Doyle, Code 6110 Washington, DC 20375
1	Commanding Officer Naval Underwater Systems Center Weapons Dept. ATTN: R.S. Lazar/Code 36301 Newport, RI 02840
2	Commander Naval Weapons Center ATTN: T. Boggs, Code 388 T. Parr, Code 3895 China Lake, CA 93555-6001
1	Superintendent Naval Postgraduate School Dept. of Aeronautics ATTN: D.W. Netzer Monterey, CA 93940
3	AL/LSCF ATTN: R. Corley R. Geisler J. Levine Edwards AFB, CA 93523-5000
1	AL/MKPB ATTN: B. Goshgarian Edwards AFB, CA 93523-5000
1	AFOSR ATTN: J.M. Tishkoff Bolling Air Force Base Washington, DC 20332
1	OSD/SDIO/IST ATTN: L. Caveny Pentagon Washington, DC 20301-7100
1	Commandant USAFAS ATTN: ATSF-TSM-CN Fort Sill, OK 73503-5600
1	F.J. Seiler ATTN: S.A. Shackelford USAF Academy, CO 80840-6528

<u>No. of Copies</u>	<u>Organization</u>
1	University of Dayton Research Institute ATTN: D. Campbell AL/PAP Edwards AFB, CA 93523
1	NASA Langley Research Center Langley Station ATTN: G.B. Northam/MS 168 Hampton, VA 23365
4	National Bureau of Standards ATTN: J. Hastie M. Jacox T. Kashiwagi H. Semerjian US Department of Commerce Washington, DC 20234
1	Aerojet Solid Propulsion Co. ATTN: P. Micheli Sacramento, GA 95813
1	Applied Combustion Technology, Inc. ATTN: A.M. Varney P.O. Box 17885 Orlando, FL 32860
2	Applied Mechanics Reviews The American Society of Mechanical Engineers ATTN: R.E. White A.B. Wenzel 345 E. 47th Street New York, NY 10017
1	Atlantic Research Corp. ATTN: M.K. King 5390 Cherokee Avenue Alexandria, VA 22314
1	Atlantic Research Corp. ATTN: R.H.W. Waesche 7511 Wellington Road Gainesville, VA 22065
1	AVCO Everett Research Laboratory Division ATTN: D. Stickler 2385 Revere Beach Parkway Everett, MA 02149

<u>No. of Copies</u>	<u>Organization</u>
1	Battelle Memorial Institute Tactical Technology Center ATTN: J. Huggins 505 King Avenue Columbus, OH 43201
1	Cohen Professional Services ATTN: N.S. Cohen 141 Channing Street Redlands, CA 92373
1	Exxon Research & Eng. Co. ATTN: A. Dean Route 22E Annandale, NJ 08801
1	Ford Aerospace and Communications Corp. DIVAD Division Div. Hq., Irvine ATTN: D. Williams Main Street & Ford Road Newport Beach, CA 92663
1	General Applied Science Laboratories, Inc. 77 Raynor Avenue Ronkonkoma, NY 11779-6649
1	General Electric Armament & Electrical Systems ATTN: M.J. Bulman Lakeside Avenue Burlington, VT 05401
1	General Electric Ordnance Systems ATTN: J. Mandzy 100 Plastics Avenue Pittsfield, MA 01203
2	General Motors Resch Labs Physics Department ATTN: T. Sloan R. Teets Warren, MI 48090
2	Hercules, Inc. Allegheny Ballistics Lab. ATTN: W.B. Walkup E.A. Yount P.O. Box 210 Rocket Center, WV 26726

<u>No. of Copies</u>	<u>Organization</u>
1	Honeywell, Inc. Government and Aerospace Products ATTN: D.E. Broden/ MS MN50-2000 600 2nd Street NE Hopkins, MN 55343
1	Honeywell, Inc. ATTN: R.E. Tompkins MN38-3300 10400 Yellow Circle Drive Minnetonka, MN 55343
1	IBM Corporation ATTN: A.C. Tam Research Division 5600 Cottle Road San Jose, CA 95193
1	IIT Research Institute ATTN: R.F. Remaly 10 West 35th Street Chicago, IL 60616
2	Director Lawrence Livermore National Laboratory ATTN: C. Westbrook M. Costantino P.O. Box 808 Livermore, CA 94550
1	Lockheed Missiles & Space Co. ATTN: George Lo 3251 Hanover Street Dept. 52-35/B204/2 Palo Alto, CA 94304
1	Los Alamos National Lab ATTN: B. Nichols T7, MS-B284 P.O. Box 1663 Los Alamos, NM 87545
1	National Science Foundation ATTN: A.B. Harvey Washington, DC 20550
1	Olin Ordnance ATTN: V. McDonald, Library P.O. Box 222 St. Marks, FL 32355-0222

<u>No. of Copies</u>	<u>Organization</u>
1	Paul Gough Associates, Inc. ATTN: P.S. Gough 1048 South Street Portsmouth, NH 03801-5423
2	Princeton Combustion Research Laboratories, Inc. ATTN: M. Summerfield N.A. Messina 475 US Highway One Monmouth Junction, NJ 08852
1	Hughes Aircraft Company ATTN: T.E. Ward 8433 Fallbrook Avenue Canoga Park, CA 91303
1	Rockwell International Corp. Rocketdyne Division ATTN: J.E. Flanagan/HB02 6633 Canoga Avenue Canoga Park, CA 91304
4	Sandia National Laboratories Division 8354 ATTN: R. Cattolica S. Johnston P. Mattern D. Stephenson Livermore, CA 94550
1	Science Applications, Inc. ATTN: R.B. Edelman 23146 Cumorah Crest Woodland Hills, CA 91364
3	SRI International ATTN: G. Smith D. Crosley D. Golden 333 Ravenswood Avenue Menlo Park, CA 94025
1	Stevens Institute of Tech. Davidson Laboratory ATTN: R. McAlevy, III Hoboken, NJ 07030
1	Thiokol Corporation Elkton Division ATTN: S.F. Palopoli P.O. Box 241 Elkton, MD 21921

<u>No. of Copies</u>	<u>Organization</u>	<u>No. of Copies</u>	<u>Organization</u>
1	Morton Thiokol, Inc. Huntsville Division ATTN: J. Deur Huntsville, AL 35807-7501	1	University of California, Berkeley Chemistry Department ATTN: C. Bradley Moore 211 Lewis Hall Berkeley, CA 94720
3	Thiokol Corporation Wasatch Division ATTN: S.J. Bennett P.O. Box 524 Brigham City, UT 84302	1	University of California, San Diego ATTN: F.A. Williams AMES, B010 La Jolla, CA 92093
1	United Technologies Research Center ATTN: A.C. Eckbreth East Hartford, CT 06108	2	University of California, Santa Barbara Quantum Institute ATTN: K. Schofield M. Steinberg Santa Barbara, CA 93106
3	United Technologies Corp. Chemical Systems Division ATTN: R.S. Brown T.D. Myers (2 copies) P.O. Box 49028 San Jose, CA 95161-9028	1	University of Colorado at Boulder Engineering Center ATTN: J. Daily Campus Box 427 Boulder, CO 80309-0427
1	Universal Propulsion Company ATTN: H.J. McSpadden Black Canyon Stage 1 Box 1140 Phoenix, AZ 85029	2	University of Southern California Dept. of Chemistry ATTN: S. Benson C. Wittig Los Angeles, CA 90007
1	Veritay Technology, Inc. ATTN: E.B. Fisher 4845 Millersport Highway P.O. Box 305 East Amherst, NY 14051-0305	1	Case Western Reserve Univ. Div. of Aerospace Sciences ATTN: J. Tien Cleveland, OH 44135
1	Brigham Young University Dept. of Chemical Engineering ATTN: M.W. Beckstead Provo, UT 84058	1	Cornell University Department of Chemistry ATTN: T.A. Cool Baker Laboratory Ithaca, NY 14853
1	California Institute of Tech. Jet Propulsion Laboratory ATTN: L. Strand/MS 512/102 4800 Oak Grove Drive Pasadena, CA 91009	1	University of Delaware ATTN: T. Brill Chemistry Department Newark, DE 19711
1	California Institute of Technology ATTN: F.E.C. Culick/ MC 301-46 204 Karman Lab. Pasadena, CA 91125	1	University of Florida Dept. of Chemistry ATTN: J. Winefordner Gainesville, FL 32611
1	University of California Los Alamos Scientific Lab. P.O. Box 1663, Mail Stop B216 Los Alamos, NM 87545		

No. of Copies	Organization
3	Georgia Institute of Technology School of Aerospace Engineering ATTN: E. Price W.C. Strahle B.T. Zinn Atlanta, GA 30332
1	University of Illinois Dept. of Mech. Eng. ATTN: H. Krier 144MEB, 1206 W. Green St. Urbana, IL 61801
1	Johns Hopkins University/APL Chemical Propulsion Information Agency ATTN: T.W. Christian Johns Hopkins Road Laurel, MD 20707
1	University of Michigan Gas Dynamics Lab Aerospace Engineering Bldg. ATTN: G.M. Faeth Ann Arbor, MI 48109-2140
1	University of Minnesota Dept. of Mechanical Engineering ATTN: E. Fletcher Minneapolis, MN 55455
3	Pennsylvania State University Applied Research Laboratory ATTN: K.K. Kuo H. Palmer M. Micci University Park, PA 16802
1	Pennsylvania State University Dept. of Mechanical Engineering ATTN: V. Yang University Park, PA 16802
1	Polytechnic Institute of NY Graduate Center ATTN: S. Lederman Route 110 Farmingdale, NY 11735

No. of Copies	Organization
2	Princeton University Forrestal Campus Library ATTN: K. Brezinsky I. Glassman P.O. Box 710 Princeton, NJ 08540
1	Purdue University School of Aeronautics and Astronautics ATTN: J.R. Osborn Grissom Hall West Lafayette, IN 47906
1	Purdue University Department of Chemistry ATTN: E. Grant West Lafayette, IN 47906
2	Purdue University School of Mechanical Engineering ATTN: N.M. Laurendeau S.N.B. Murthy TSPC Chaffee Hall West Lafayette, IN 47906
1	Rensselaer Polytechnic Inst. Dept. of Chemical Engineering ATTN: A. Fontijn Troy, NY 12181
1	Stanford University Dept. of Mechanical Engineering ATTN: R. Hanson Stanford, CA 94305
1	University of Texas Dept. of Chemistry ATTN: W. Gardiner Austin, TX 78712
1	University of Utah Dept. of Chemical Engineering ATTN: G. Flandro Salt Lake City, UT 84112
1	Virginia Polytechnic Institute and State University ATTN: J.A. Schetz Blacksburg, VA 24061

<u>No. of</u> <u>Copies</u>	<u>Organization</u>
1	Freedman Associates ATTN: E. Freedman 2411 Diana Road Baltimore, MD 21209-1525

INTENTIONALLY LEFT BLANK.

USER EVALUATION SHEET/CHANGE OF ADDRESS

This Laboratory undertakes a continuing effort to improve the quality of the reports it publishes. Your comments/answers to the items/questions below will aid us in our efforts.

1. BRL Report Number BRL-TR-3163 Date of Report OCTOBER 1990
2. Date Report Received _____
3. Does this report satisfy a need? (Comment on purpose, related project, or other area of interest for which the report will be used.) _____

4. Specifically, how is the report being used? (Information source, design data, procedure, source of ideas, etc.) _____

5. Has the information in this report led to any quantitative savings as far as man-hours or dollars saved, operating costs avoided, or efficiencies achieved, etc? If so, please elaborate. _____

6. General Comments. What do you think should be changed to improve future reports? (Indicate changes to organization, technical content, format, etc.) _____

CURRENT
ADDRESS

Name

Organization

Address

City, State, Zip Code

/ If indicating a Change of Address or Address Correction, please provide the New or Correct Address in Block 6 above and the Old or Incorrect address below.

OLD
ADDRESS

Name

Organization

Address

City, State, Zip Code

(Remove this sheet, fold as indicated, staple or tape closed, and mail.)

-----FOLD HERE-----

DEPARTMENT OF THE ARMY

Director
U.S. Army Ballistic Research Laboratory
ATTN: SLCBR-DD-T
Aberdeen Proving Ground, MD 21005-5066
OFFICIAL BUSINESS




**NO POSTAGE
NECESSARY
IF MAILED
IN THE
UNITED STATES**

BUSINESS REPLY MAIL
FIRST CLASS PERMIT No 0001, APG, MD

POSTAGE WILL BE PAID BY ADDRESSEE

Director
U.S. Army Ballistic Research Laboratory
ATTN: SLCBR-DD-T
Aberdeen Proving Ground, MD 21005-9989

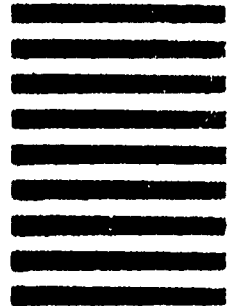

-----FOLD HERE-----